




c) Choose the destination folder and then press “Next”;

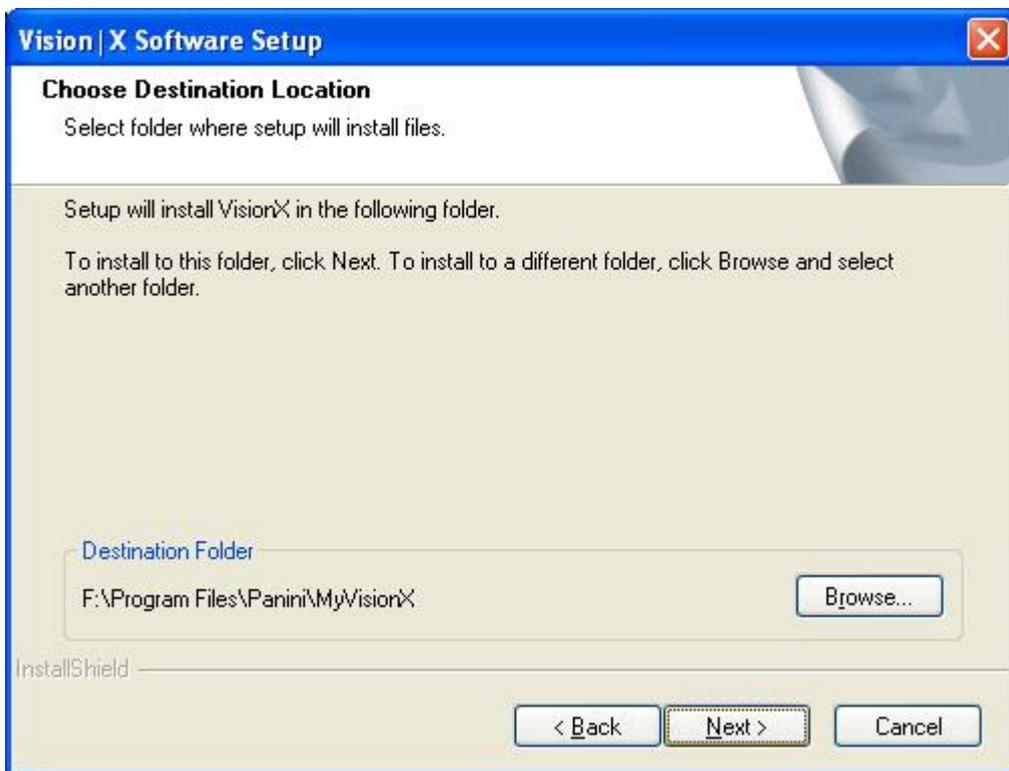

d) Select the desired features and then press “Next”;

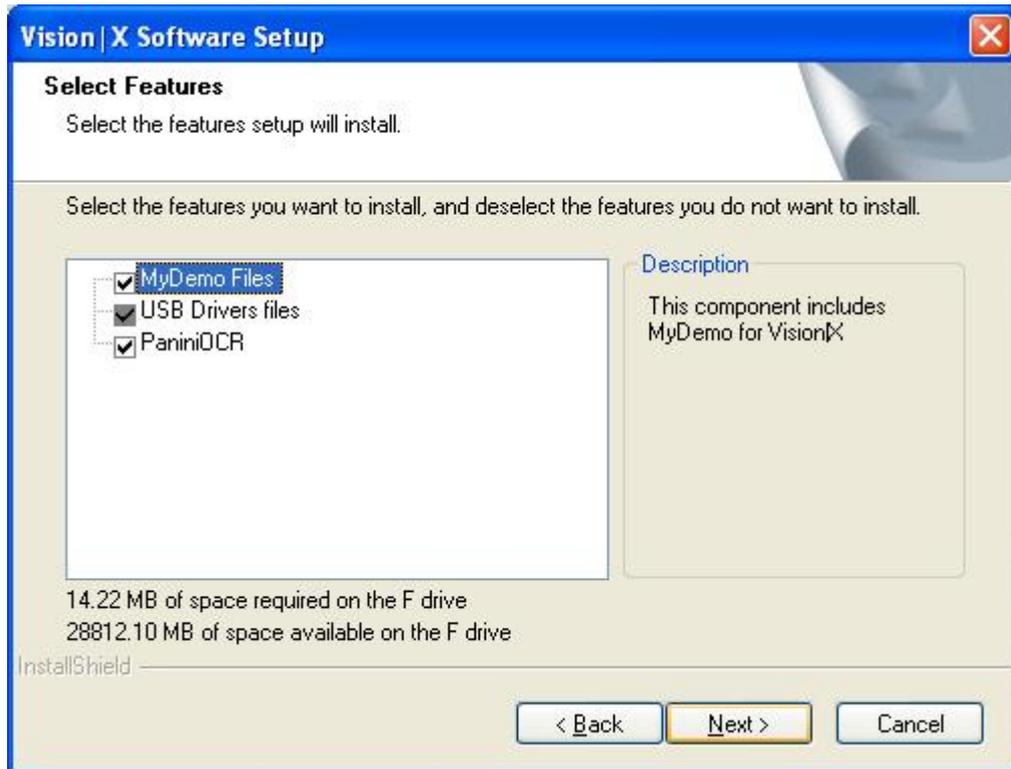

e) Press "Finish" to complete the setup.

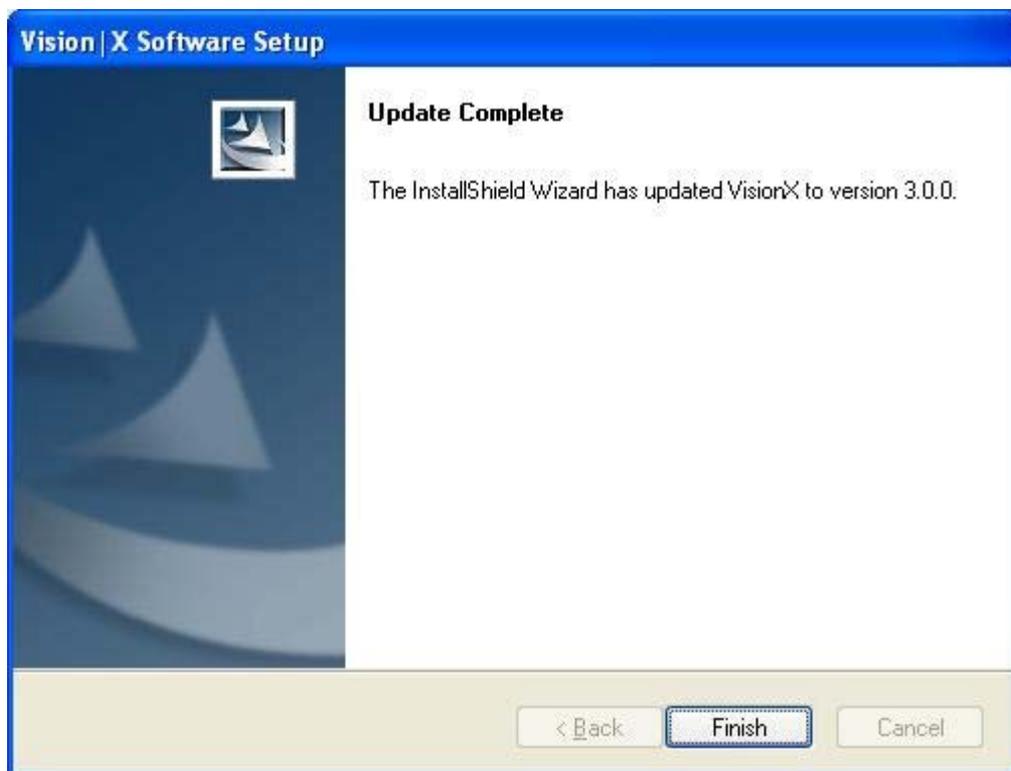

It's also possible to make these operations in "silent" mode. In order to enable this option the user has to launch the setup.exe application using the following command line: "setup silent /hide_splash /hide_progress".



This command runs the upgrade avoiding the user intervention, using all the default choices.

It's also possible to force the remove of a previous version. In order to remove a previous version use the following command line:

“setup /removeonly”.

3. Uninstall (the same Vision X software version installed).

The 1st time the user runs the setup, it requests him to remove the installed version.

The 2nd time the user runs the setup it installs the software.